

Configurable Middleware



Integrate
in-vitro diagnostic
Instruments with LIS/EHR/HIS

INDUSTRIES

- Referral Labs
- Hospital Labs
- Clinical Research
- Physician's Office Lab
- Bio Banking
- Blood Banks
- Clinical Research
- Molecular Diagnostics



Interfacer Features

- ✓ Middleware for instrument interfacing.
- ✓ Middleware that simplifies communication between any LIS/EHR/HIS & clinical diagnostic Instruments.
- ✓ Fully configurable communication protocol (ASTM 1381, ASTM 1394, HL7 and custom protocols).
- ✓ Supports both Uni-directional & Bi-directional communication with instrument.
- ✓ Real-time response and full instrument monitor and control.
- ✓ Advanced integration with any LIS/EHR/HIS through HL7 and other methods.
- ✓ One time mapping of Tests and Test parameters between LIS/EHR/HIS & Instruments.
- ✓ H IPAA compliant.
- ✓ Secure communication.

Simple & User Friendly Design

Interfacer is a simple and user friendly middleware designed for interfacing in-vitro diagnostic clinical instruments in hospitals, referral laboratories and research facilities.


Interfacer acts a middleware between your LIS/EHR/HIS and the clinical laboratory instruments, Runs as a service in server and has inbuilt features to connect and configure any clinical instrument irrespective of the protocol via both TCP/ IP and RS232 modes

Driverless Technology

Traditional middleware solutions adopted drivers for specific manufacturer, instrument model and software version. Interfacer is a driverless solution. Interfacer can be setup and configured to communicate with any IVD instrument from any manufacture, model or software version.


The communication with instrument is fully transparent and the system configurator can easily configure the interface with the instrument.

- 1 Service oriented middleware architecture for easy deployment, with fully scalable and flexible enterprise class middleware designed to streamline high throughput instrument interfacing.
- 2 Easy Integration of LIS/EHR/HIS or other healthcare systems in a seamless manner.

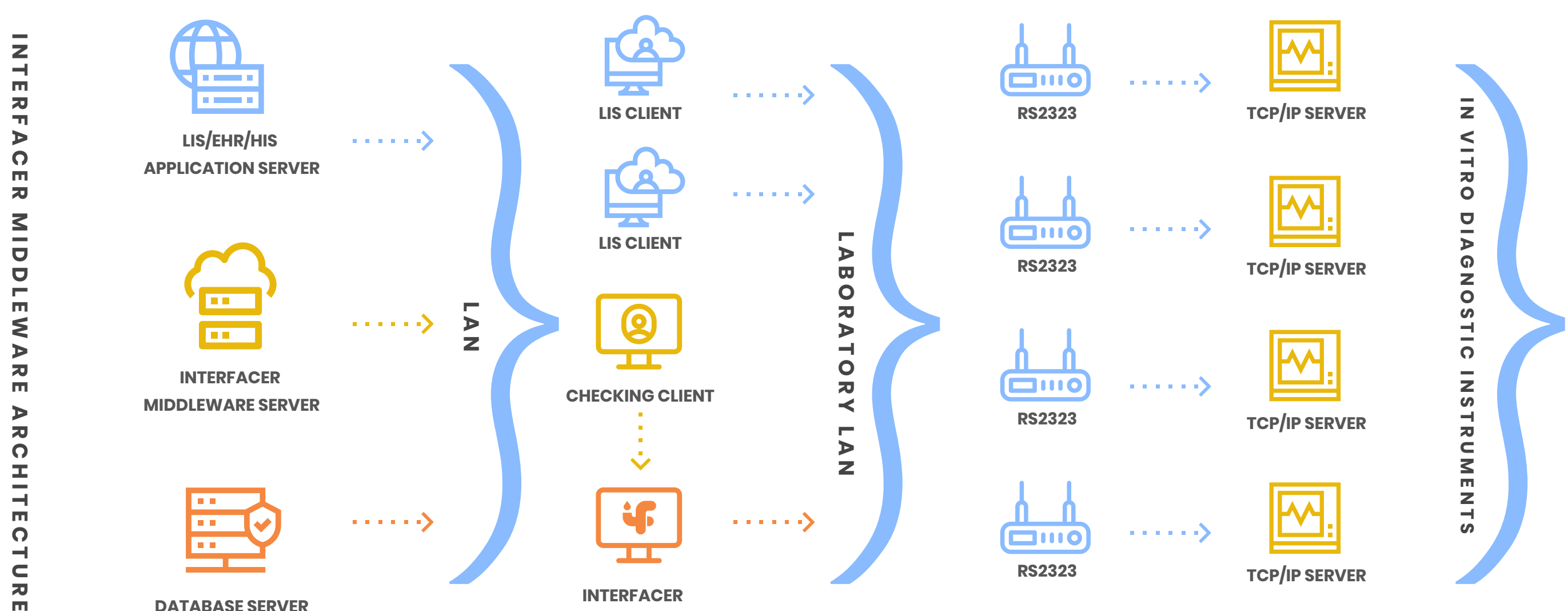

Interfacer

Bi-Directional Communication

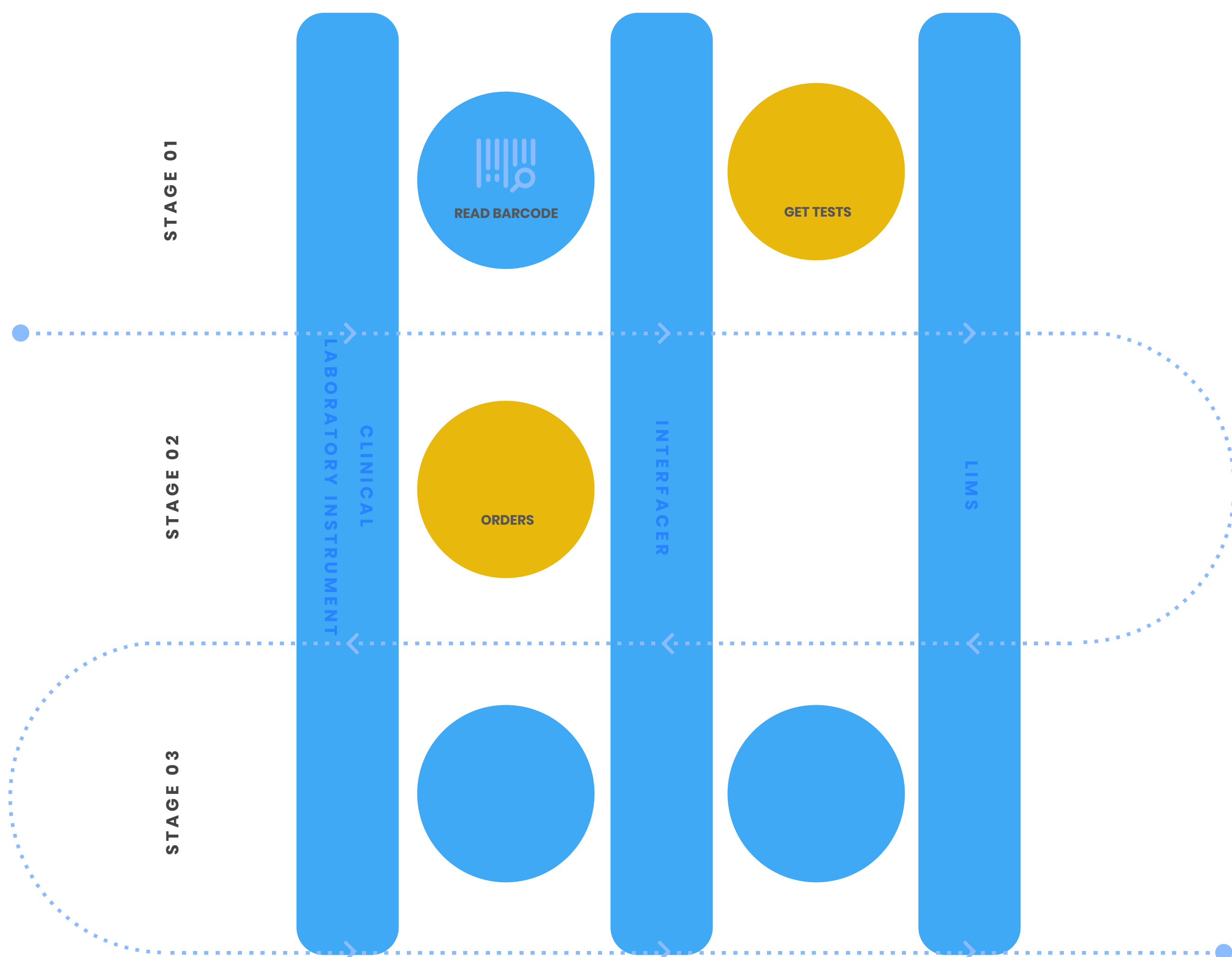

Uni-Directional & Bi-Directional Interface


Interface receives results from instruments and extracts individual tests and Parameter results and posts the results to a middleware database table, creates files (CSV, txt, xml) or HL7 messages for the LIS/EHR/HIS to pickup.

Interface communicates with the LIS/EHR/HIS through database, files or HL7. Orders placed in LIS/EHR/HIS is received in the middleware for all instruments connected to the middleware.

Instrument when loaded with the sample with bar-code enquires the middleware for test orders. middleware places the order for tests for samples. The instrument will return the results to the middleware once the results are ready. The whole communication happens seamless. The middleware/LIS mapped results are sent back to the LIS/EHR/HIS.

Technical Specification

Feature	Details
Instrument Supported	Any In vitro diagnostic instrument used in haematology, biochemistry, immunodeficiency-assay microbiology, histopathology, POCT, urine, blood gas, electrochemical, molecular diagnostic, sample procession robotics.
Manufactures supported	Any manufacturer who supports input/output via RS232, TCP/IP, USB, File, HL7.
Communication protocols supported (instrument)	ASTM 1394, ASTM 1381, HL7, Custom protocols.
Database used by middleware	MS-SQL Server (2008 or higher).
LIS/EHR/HIS Integration	Database, File based (CSV, xml, txt), HL7.
Database supported – LIS	MS-SQL Server (2008 or higher), Oracle 10g or higher, PostgreSQL 8 or higher.
Database support – middleware side	MS-SQL Server (Express, Standard, enterprise editions).
Compliance	21 CFR Part 11, HIPAA, cGCP, GLP.
Installation and Support availability	Global.
Middleware Instrument configuration	Instrument setup, communication protocol, multi-messages per instrument, message parsing, parsed field mapping
Middleware LIS/EHR/HIS configuration	Test Mapping, Parameter Mapping, Results Mapping, Result flagging.
EDI	Push and Pull of files for communication.

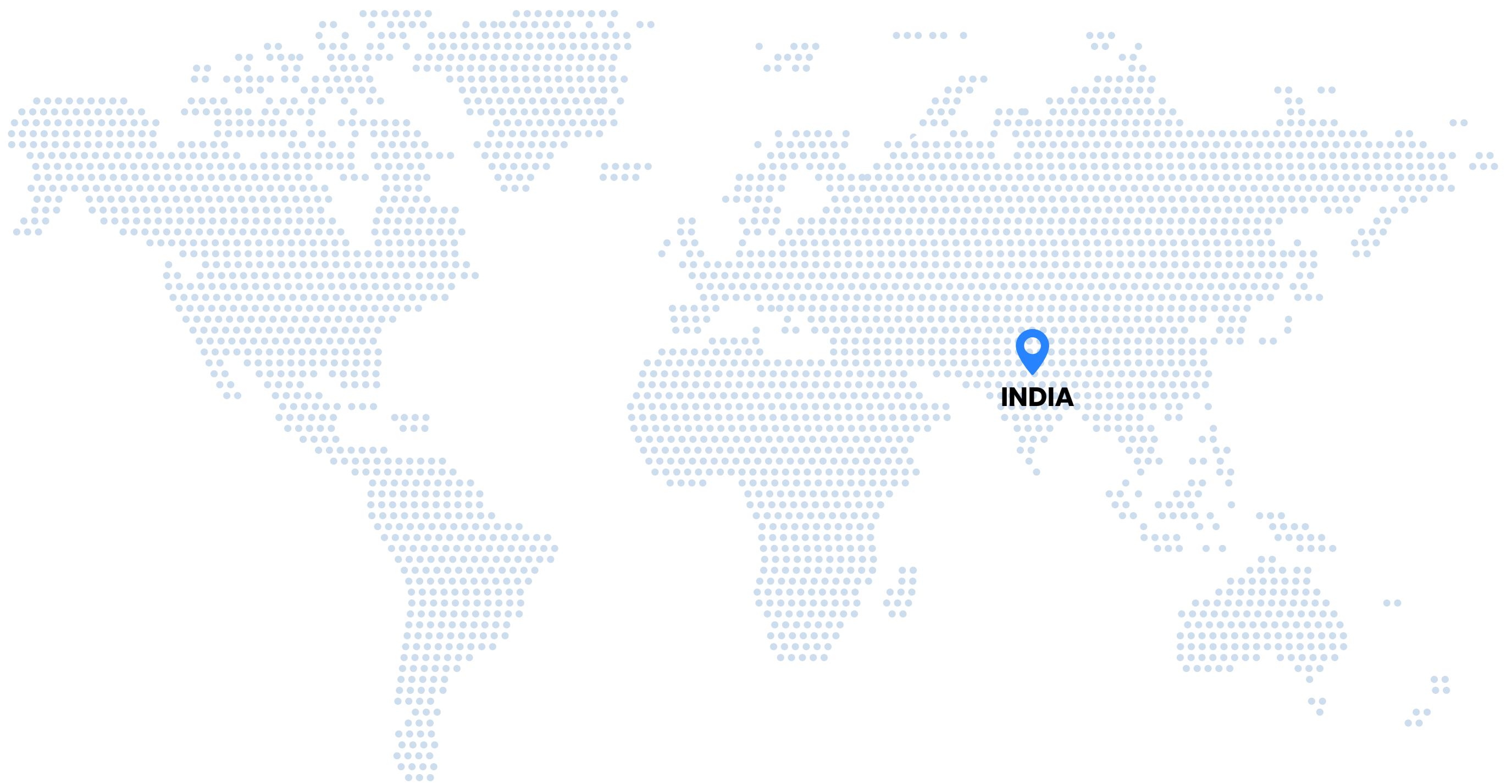

Locations

Global HQ

Priyan Plaza, 76, Nelson Road,
Aminjikarai, Chennai 600 029, INDIA
T: +91 44 4208 2005
T: +91 44 42189406
E: info@agaramtech.com

Our Products



LIMS

Qualis



CFR

Gateway



ELN

Logilab



SDMS

Logilab



Interfacer



STABILITY

Qualis

Disclaimer:

www.agaramtech.com ©2020 Agaram Technologies. All rights reserved. For information regarding this and other Agaram Technologies products, please contact Agaram Technologies.

The information contained in this document is proprietary and confidential to Agaram Technologies. No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, including photocopying and recording, for any purpose, without the express written permission of Agaram Technologies.

This document is subject to change without notice, and Agaram Technologies does not warrant that the material contained in this document is error-free. If you find any problems with this document, please report them to Agaram Technologies in writing. Qualis and Logilab are trademarks of Agaram Technologies. All other trademarks are property of their respective owners.

This document may contain statements concerning possible functionality for Qualis and Logilab software products. Agaram Technologies disclaims any express or implied commitment to deliver functionality or software unless or until generally commercially available. Any statements of possible future direction are for information purpose only and Agaram Technologies makes no express or implied commitments or representations concerning the timing and content of any future functionalities or releases.